


# Sialendoscopy

SALIVARY DUCT ACCESS & INTERVENTION

*Minimally invasive solutions*  
designed for gland preservation

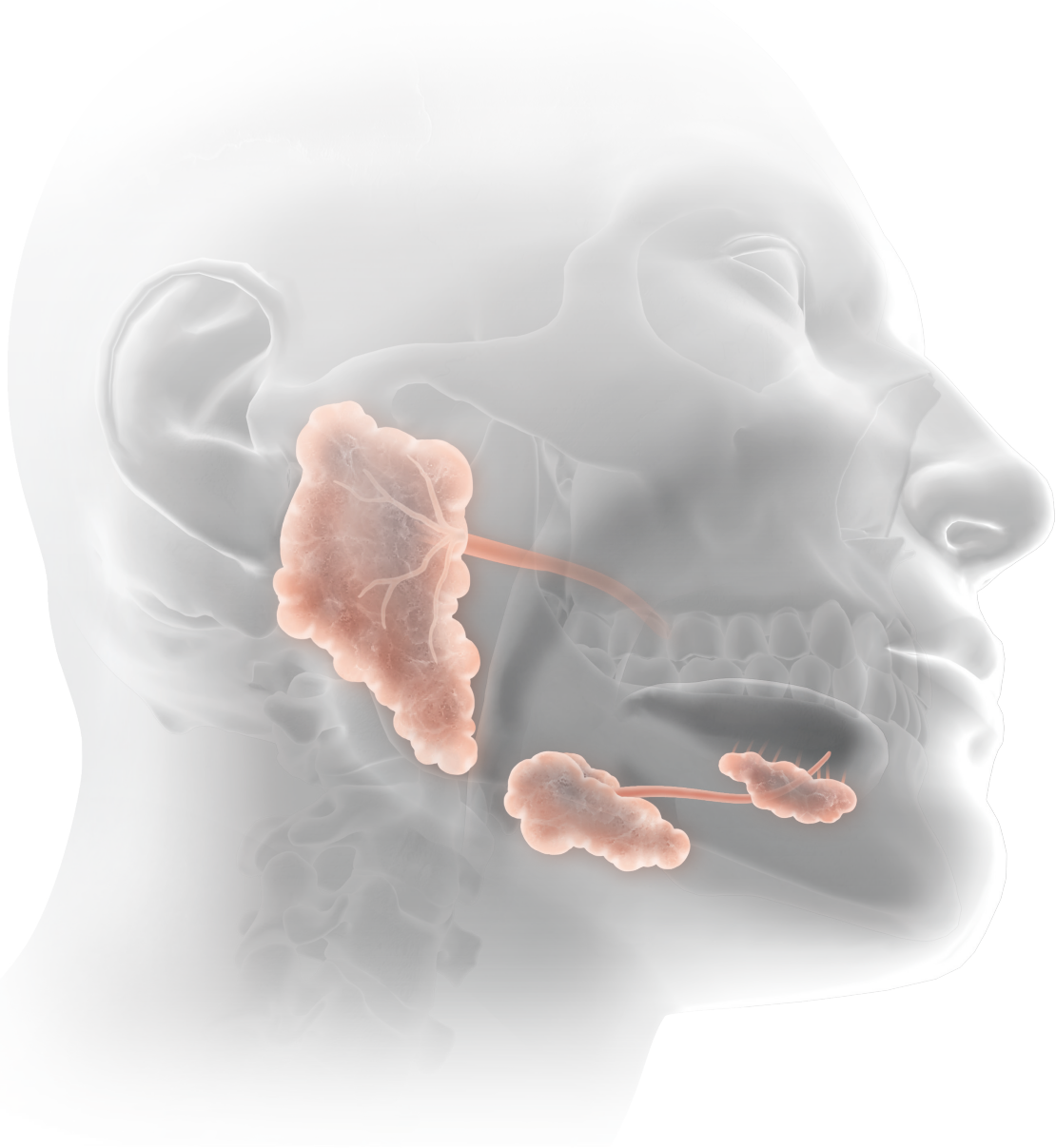

Access & Treatment Essentials

Everis™


# Interventional Sialendoscopy

The minimally invasive treatment of obstructive salivary gland disease

ACCESS  
PRODUCTS →

INTERVENTION  
PRODUCTS →

ORDERING  
INFORMATION →

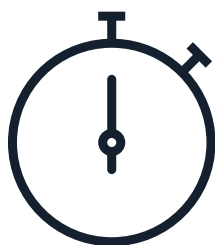

REDUCED  
RECOVERY TIME



NO  
INCISION



GLAND  
PRESERVATION

Sialendoscopy offers a variety of advantages to the patient, including reduced recovery time, no incision or subsequent scarring, minimized potential for nerve damage, preservation of the salivary gland, and early return to a normal diet.<sup>3, 12</sup>

Sialendoscopy allows stones and strictures to be managed in office with a low risk of major complications to the patient compared to open surgical procedures.<sup>1</sup>

\*A combined approach may be necessary depending on stone size, mobility, and clinical judgment.

Salivary duct access tools allow me to gain and maintain access to the duct, which gives me much greater control throughout the entire procedure.

**Jack Kolenda, MD**


Chief, Department of Otolaryngology  
Oakville-Trafalgar Memorial Hospital  
Oakville, Ontario

Everis continually collaborates with physicians to improve procedures, including partnering with leading physicians to create a comprehensive line of salivary duct access and interventional products that are designed to maintain access throughout the procedure and treat obstructive salivary gland disorders. This line of products is designed to help reduce trauma for the patient and increase efficiency with every step.

*Dr. Kolenda is the inventor of the Kolenda Salivary Access Introducer Set and received a royalty payment to use this technology.*


## ACCESS

**Gain and maintain access to the salivary duct throughout the entire procedure.**

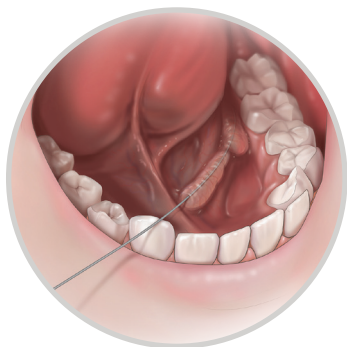


### Wire Guide

Access the salivary duct through the papilla.

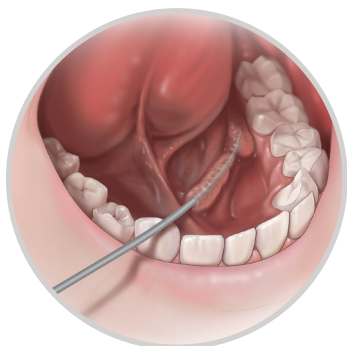


### Salivary Access Dilator Set

Introduce a series of flexible dilators over the wire guide to expand the papilla and prepare the salivary duct for the introduction of procedural instruments.

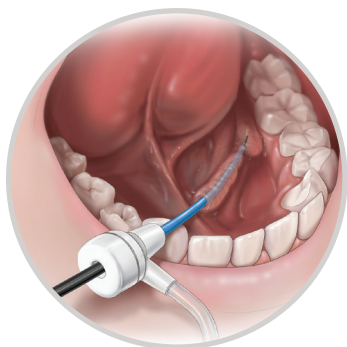


### Kolenda Salivary Access Introducer Set

Create an open working channel into the salivary duct by passing the sheath over the wire guide. The sheath is also designed to protect the ductal wall as you insert and remove procedural instruments.

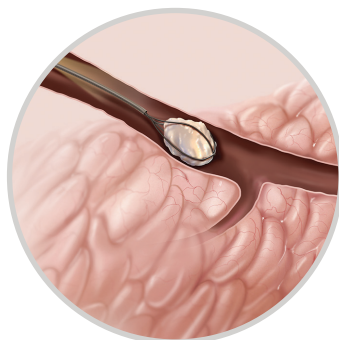


### NCircle® Tipless Salivary Stone Extractor

Insert the NCircle® nitinol basket into the duct to remove salivary stones and fragments. The stone extractor has a tipless configuration that allows you to position the basket against the mucosal lining.

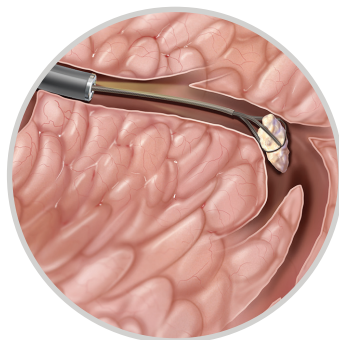


### NGage® Salivary Stone Extractor

Insert the NGage® nitinol basket into the duct to remove salivary stones and fragments. The stone extractor works as both a basket and a grasper.

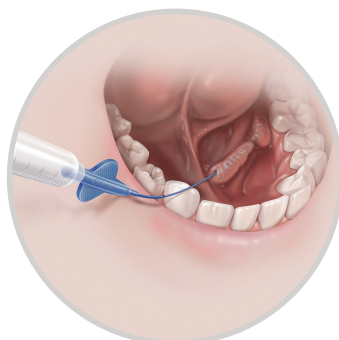

### SialoCath® Salivary Duct Catheter

Insert the catheter over a wire guide or through the Kolenda Salivary Access Introducer to irrigate the submandibular or parotid salivary ducts to flush stone fragments.

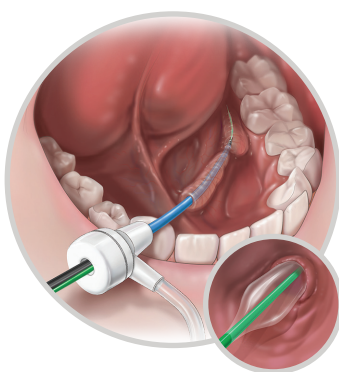

### Advance® Salivary Duct Balloon Catheter Set

Once you have created a working channel in the duct, insert the balloon catheter over a wire guide and use it alongside a sialendoscope for direct visualization to dilate at the stricture site.

## ORDERING INFORMATION

### Wire Guide


| Order Number | Reference Part Number | Wire Guide Diameter<br>inch | Wire Guide Length<br>cm | Tip Configuration |
|--------------|-----------------------|-----------------------------|-------------------------|-------------------|
| G02274       | C-SF-15-50            | 0.015                       | 50                      | Straight          |

### Salivary Access Dilator Set

| Order Number | Reference Part Number | Dilator Fr | Dilator OD<br>mm | Dilator OD<br>inch |
|--------------|-----------------------|------------|------------------|--------------------|
| G24131       | SDDS-18-20-HC         | 4.0        | 1.40             | 0.055              |
|              |                       | 5.0        | 1.67             | 0.066              |
|              |                       | 6.0        | 2.0              | 0.079              |
|              |                       | 7.0        | 2.33             | 0.092              |

### Kolenda Salivary Access Introducer Set

| Order Number | Reference Part Number | Introducer ID<br>mm/inch | Wire Guide Diameter<br>inch | Sheath Length<br>cm |
|--------------|-----------------------|--------------------------|-----------------------------|---------------------|
| G24129       | SDAS-6.0-18           | 2.21/0.087               | 0.018                       | 5.5                 |
| G24132       | SDAS-5.0-18           | 1.88/0.074               | 0.018                       | 5.5                 |

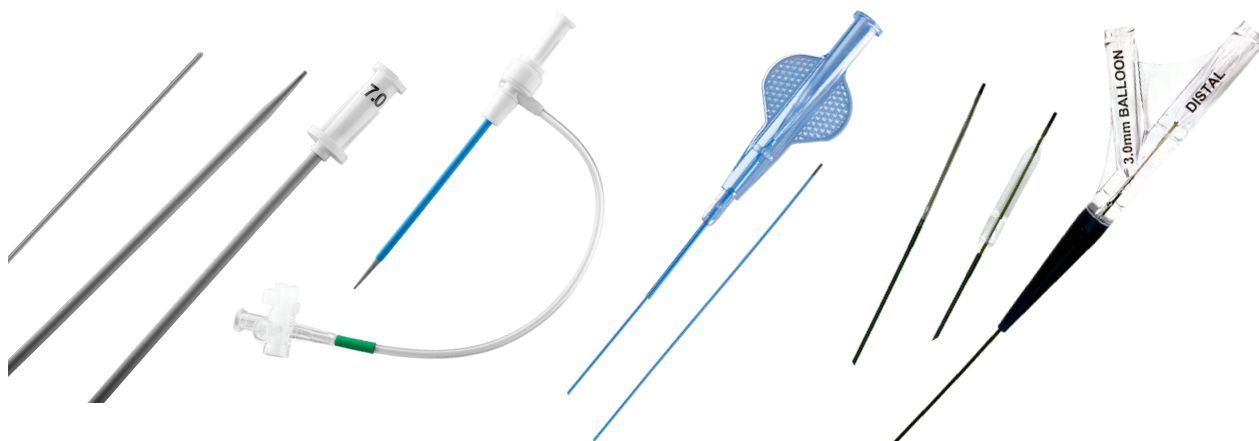

## NGage® Salivary Stone Extractor

| Order Number | Reference Part Number | Catheter Fr | Catheter Length cm | Basket Diameter mm |
|--------------|-----------------------|-------------|--------------------|--------------------|
| G23798       | SSEB-1.7-115-8        | 1.7         | 115                | 8                  |

## NCircle® Tipless Salivary Stone Extractor

| Order Number | Reference Part Number | Catheter Fr | Catheter Length cm | Basket Diameter mm |
|--------------|-----------------------|-------------|--------------------|--------------------|
| G23797       | SSEB-1.5-115-10       | 1.5         | 115                | 10                 |

## SialoCath® Salivary Duct Catheter

| Order Number | Reference Part Number | Catheter OD mm/inch | Catheter ID mm/inch | Compatible Wire Guide Diameter inch | Compatible Catheter Length cm | Tip Configuration |
|--------------|-----------------------|---------------------|---------------------|-------------------------------------|-------------------------------|-------------------|
| G35166       | SDC-2.2-20            | 0.74/0.029          | 0.53/0.021          | 0.018                               | 15                            | Straight          |

## Advance® Salivary Duct Balloon Catheter Set

| Order Number | Reference Part Number | Wire Guide Diameter inch | Inflated Balloon Diameter mm |
|--------------|-----------------------|--------------------------|------------------------------|
| G34991       | SDBS-14-30-1.5-2      | 0.014                    | 1.5                          |
| G34992       | SDBS-14-30-2.0-2      | 0.014                    | 2.0                          |
| G34993       | SDBS-14-30-2.5-2      | 0.014                    | 2.5                          |
| G34994       | SDBS-14-30-3.0-2      | 0.014                    | 3.0                          |

## References

1. Nahlieli O, Nakar LH, Nazarian Y, et al. Sialendoscopy: a new approach to salivary gland obstructive pathology. *J Am Dent Assoc*. 2006;137(10):1394–1400.
2. Carta F, Farneti P, Cantore S, et al. Sialendoscopy for salivary stones: principles, technical skills and therapeutic experience. *Acta Otorhinolaryngol Ital*. 2017;37(2):102–112.
3. Achim V, Light TJ, Andersen PE. Gland preservation in patients undergoing sialendoscopy. *Otolaryngol Head Neck Surg*. 2017;157(1):53–57.

® indicates U.S. trademark registration. All trademarks and/or images are the property of their respective owners or holders.



## Sialendoscopy

SALIVARY DUCT ACCESS & INTERVENTION

For more detailed product information, including the Instructions for Use, which contains the indication statement, contraindications, precautions, and potential complications, visit [www.everismedical.com](http://www.everismedical.com).



### Everis

555 East Eliza St, Ste. A  
Schoolcraft, MI 49087 USA  
t: 888 902 2239  
[www.everismedical.com](http://www.everismedical.com)

DOC-6570 Rev 1  
Copyright © 2026 Everis

# Everis™