
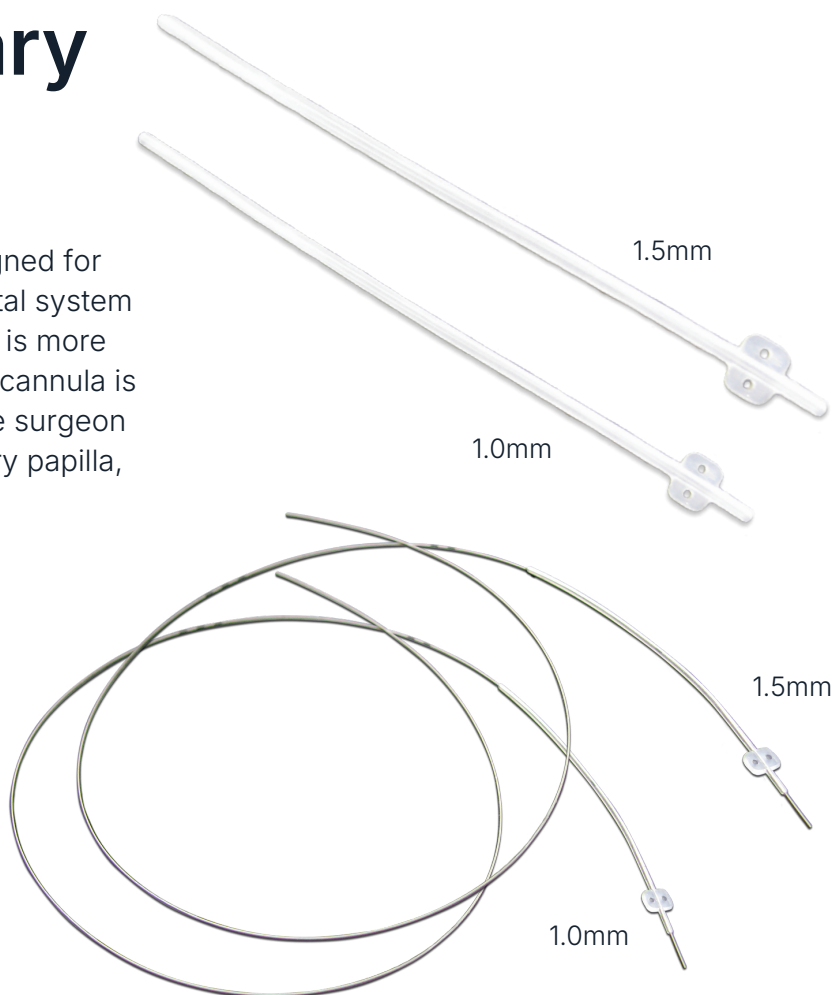


Schaitkin Salivary Duct Cannula


The Schaitkin Salivary Duct Cannula is designed for the short-term intubation of the salivary ductal system and to hold open the ductal tissue. The tube is more than sufficient length for most patients. The cannula is placed at the end of the procedure when the surgeon is concerned about the healing of the salivary papilla, Wharton's duct, or Stenson's duct.

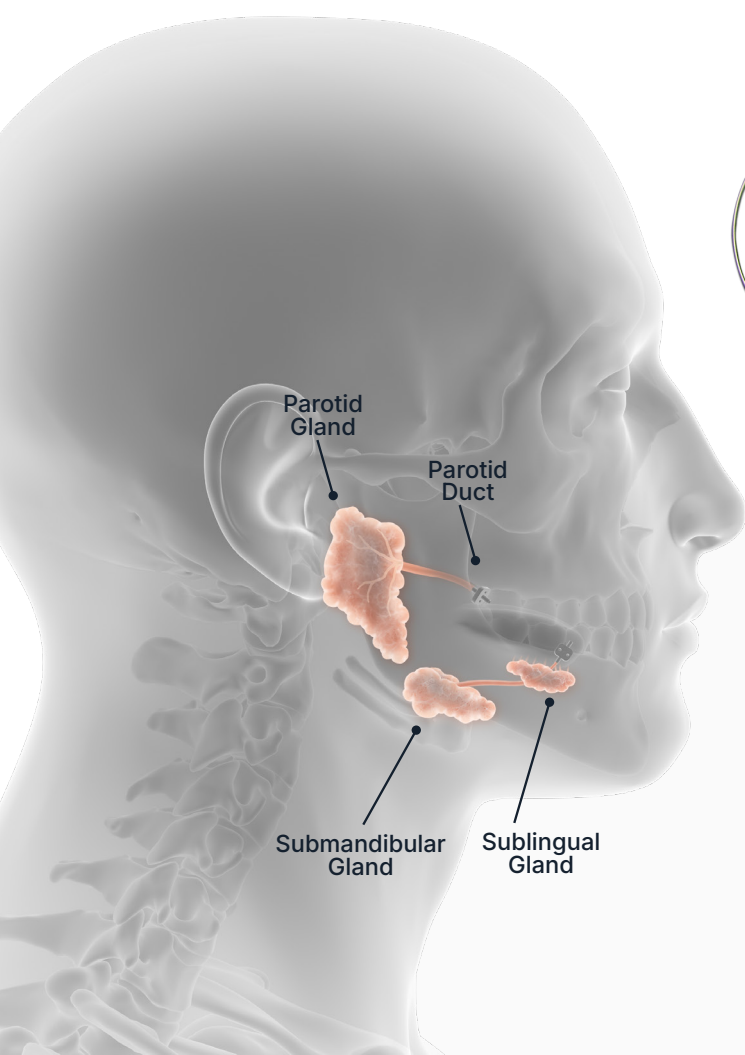


Lumen allows the use of the **Guide Wire** to ease introduction into the papilla, Wharton's duct, or Stenson's duct.

Benefits

- **Short-term stenting** of the salivary papilla, Wharton's duct, or Stenson's duct.
- **Reduces chances of scarring** post endoscopic stone removal.
- Allows **continued saliva flow** through ducts.
- **Convenient suture points** for securing after placement.
- Made of **soft implant-grade silicone** that is well tolerated by patients.

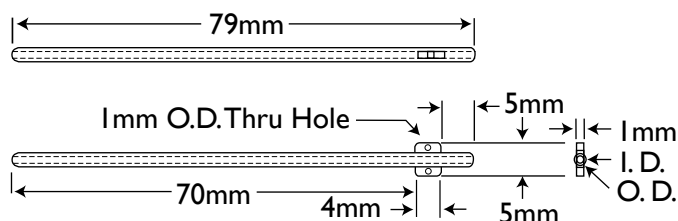

Schaitkin Salivary Duct Cannula

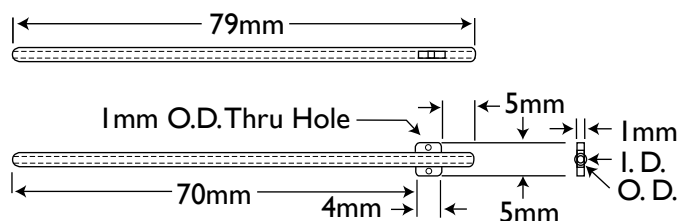

Schaitkin Salivary Duct Cannula with Guide Wire

Ordering Information

Part Number	Outside Diameter	Inside Diameter	Length
Schaitkin Salivary Duct Cannula			
SDC-2.0-S-ST (STERILE)	2.0mm	0.78mm	
Guide Wire			
GW-018-ST	0.45mm		45cm
Schaitkin Salivary Duct Cannula with Guide Wire			
SDC-1.0-S-GW	1.0mm	0.5mm	
SDC-1.5-S-GW	1.5mm	0.53mm	

Learn more at www.everismedical.com

To place an order,
email sales@everismedical.com

Ultra-smooth^{Plus}™

All products have been treated with Ultra-smooth Plus[®] surface treatment. This proprietary technology modifies the surface properties of silicone. Ultra-smooth Plus[®] treated silicone is thromboresistant, resistant to biofilm germination, bacterial and fungal growth, and has less surface friction. Advantages of Ultra-smooth Plus[®] include:

- Slippery surface—smooth, less friction than untreated products, no need for lubricants
- Biocompatibility—thromboresistant, resistant to biofilm formation, resistant to bacterial and fungal growth
- Extended product life—improves hydrophilicity, wear resistant, abrasion resistant
- Reduced rates of infection and thrombosis

INDICATIONS: Patients who have had more than basic dilation of the papilla may benefit from short-term intubation to allow for drainage of the salivary duct. The cannula can be placed via a sialodochotomy. Patients who have had parotid procedures done externally may have a retrograde insertion.

Please refer to the Intended Use and Instruction Manual for detailed information on insertion technique, sterilization, and postoperative care.

© indicates U.S. trademark registration. All trademarks and/or images are the property of their respective owners or holders.

Everis

555 East Eliza St, Ste. A
Schoolcraft, MI 49087 USA
t: 888 902 2239
www.everismedical.com

DOC-14367 Rev 3
Copyright © 2026 Everis